User Manual

Bathroom Scale

Please read this User Manual carefully before using the device

Factory: SHENZHEN SHINE INDUSTRIAL CO., LTD.
Address: 2&3F, Bldg5, 1st Industrial Zone, Changzhen Community,
Yutang Street, Guangming District, Shenzhen City, Guangdong
Province, China, 518132
E-mail: sales03@iascent.com.cn

Company: Apex CE Specialists GmbH
Address: Grafenberger Allee 277, 40237 Düsseldorf, Westfalen, Deutschland
Zip Code: 40237
Contact: Wells Yan
E-mail: info@apex-ce.com
Tel: +491776995991

Company: APEX CE SPECIALISTS LIMITED
Address: 6f East Union Street, Manchester, England, M16 9AE
Contact: Wells
E-mail: Info@apex-ce.com
Tel: +441616371080

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09-10 Italiano


First Time Setup

IMPORTANT! BEFORE USE: Remove the isolator tab from the battery compartment


Weighing Operation


- For best accuracy, place the scale on a hard and flat floor. Carpeted or uneven floors can negatively effect the accuracy .
- The scale operates with step-on technology but will be calibrated automatically under following occasions.
A. when first time use, moving the scale or battery exchange, the scale will display scale itself weight, and appear "CAL" after around 1second, and then display "0.0" after around 2seconds , now it's ready to step on to weigh normally.

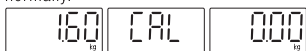B. If step on the scale before calibration, the scale will display an wrong weight and appear "CAL" after around 1 second, and then display "0.0" after around 2seconds.

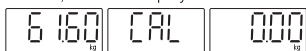
After calibration, please step on again to weigh and get correct weight.

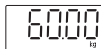
- Under other occasions, please step on to weigh directly.
- When weighing, please remain standing evenly on the scale without shifting feet or body to get right weight.
- Step off the scale, the weight readings remain for several seconds. Then the scale will turn off automatically.

VERY IMPORTANT NOTICE: If the scale placed in a **NON-HORIZONTAL POSITION OR MOVED BEFORE WEIGHING**, please note there are two situations of measurements after stepping off the scale: 1.If "CAL" displayed after around 1 second and then display "0.0" after around 2 seconds, it means that the displayed weight is incorrect, please step on the scale again to get the correct weight; 2. If there is no "CAL" display after around 2 seconds, the displayed weight is correct.

Error Messages

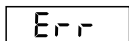If "Err" is displayed, the scale would be overloaded.
If "Err" is displayed, the weighing procedure failed, please repeat.

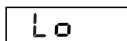When you step on the scale to turn it on, there is no display at all or "Lo" is displayed, try new batteries.

Safety Cautions

Keep the scale off water

Don't put it on uneven floor while weighing

Use a wet cloth To clean the scale

Keep the scale with care

Don't step on the edge of the scale

Slippery when wet

Specification

Capacity: Min. 5 kg/11lb/0.11lb
Max.182kg/400lb/28st

Weight unit: kg/lb/st

Graduation: 0.05kg/0.1lb/0.1lb

Dimensions: 280x280x24mm or 277x277x24mm(Depended on the Mode)

Platform: Tempered glass

Display: LCD or LED (Depended on the model)

Batteries: 2xAAA or 3xAAA (Depended on the model)

Warranty Information

VitaFit warrants all products to be of the highest quality in material, craftsmanship and service for a warranty of 3 years, effective from the date of purchase.

For more information please visit www.51vitaFit.com, www.vitaFit.com

Email: service@51vitaFit.com

Warning

1. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
2. Children shall not play with the appliance.
3. Cleaning and user maintenance shall not be made by children without supervision.
4. Model of battery that appliance used is 2xAAA or 3xAAA.
5. Non-rechargeable batteries are not to be recharged.
6. Different types of batteries or new and used batteries are not to be mixed.
7. Batteries are to be inserted with the correct polarity.
8. Exhausted batteries are to be removed from the appliance and safely disposed of.
9. If the appliance is to be stored unused for a long period, the batteries are removed.
10. The supply terminals are not to be short-circuited.
11. Take the scale, turn it over, open the battery compartment, put 2xAAA or 3xAAA batteries according to the right polarities, then return the battery compartment cover.
12. Batteries, when swallowed, can be life-threatening. You should therefore keep the batteries and scale out of reach of small children. If a battery has been swallowed, medical assistance must be sought immediately.
13. Keep packaging material away from children. There is a risk of choking if swallowed. Also pay attention that children do not pull the packaging material over their heads.
14. Avoid contact of the battery fluid with skin, eyes, and mucosa. In case of contact, immediately rinse the affected areas with plenty of clean water and consult a doctor.

Disposing of the product

Dispose of the product in accordance with the regulations that apply in your country.

Devices must not be disposed of in normal household waste.

This marking indicates that this product should not be disposed of other household waste throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmentally safe recycling.

Batteries and rechargeable batteries should not be disposed of in household waste.

All batteries and rechargeable batteries must be brought to a local or district collection point or to the retailer. Such batteries and rechargeable batteries can thus be disposed of in an environmentally friendly way.

Shenzhen Shine Industrial Co., Ltd. hereby declares that this (Name: Body scale, Model: VT725, VT7251, VT7252, VT7253, VT7255, VT7256, VT7257, VT7258, VT7259, VT726, VT7261, VT7262, VT7263, VT7265, VT7266, VT7267, VT7268, VT7269, VT727, VT7271, VT7272, VT7273, VT7275, VT7276, VT7277, VT7278, VT7279, VT1703, VT17031, VT17032, VT17033, VT17035, VT17036, VT17037, VT17038, VT17039) is compliance with the essential requirements and other relevant provisions of Directive 2014/30/EU. This declaration will lose its validity if any modification is made to the scale without our consent.

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